

ELLISYS SENSE

Fractional Microneedle RF System

ACNE




Crow feet,
Eyebag



Wrinkle



Body contouring /
Lifting

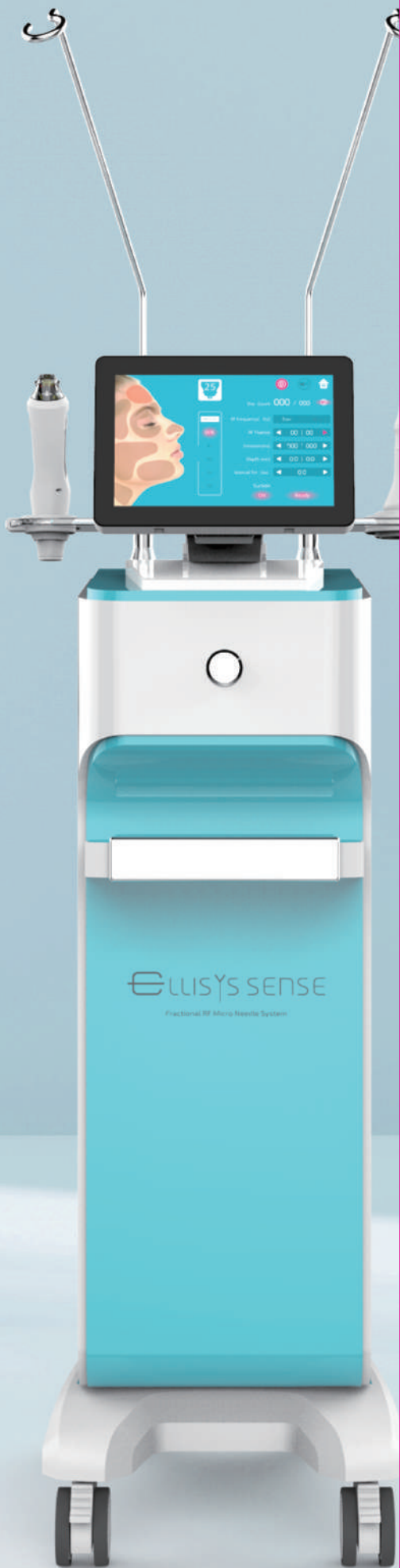

ELLISYS SENSE

Fractional Microneedle RF System

FDA

Specification

Model	ELLISYS SENSE (CWM-930S)
RF Level	0~10
RF Frequency	1Mhz & 2 Mhz
Length	Max 5mm for Mono-polar / Max 3.5mm for Bi-polar
RF output	10~900ms
Dimension	365(W) x 505(D) x 1555(H) mm

ELLISYS SENSE Remove Unwanted Skin imperfections

ELLISYS SENSE is a Fractional Microneedle Radiofrequency device that excels in tissue coagulation using electronically controlled microneedles. By delivering RF energy to the epidermis, dermis, and subcutaneous layers, it effectively heats targeted skin and controls electrothermal effects. It is particularly effective for treating scars, enlarged pores, wrinkles, uneven skin tone, saggy skin, acne, striae distensae, syringoma, milia, blackheads, and lower eyelid fat.


Monopolar RF

sends RF energy through the body to a grounding pad, heating tissues deeply.

Bipolar RF


uses two close electrodes on the device to heat specific tissues between them.

Dual-Depth Double Shot

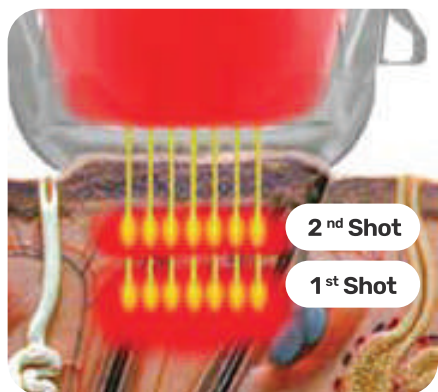


RF energy emits twice from different layers of skin in a single shot. Experience excellent skincare with one pass while reducing pain and saving time.

Vacuum Suction

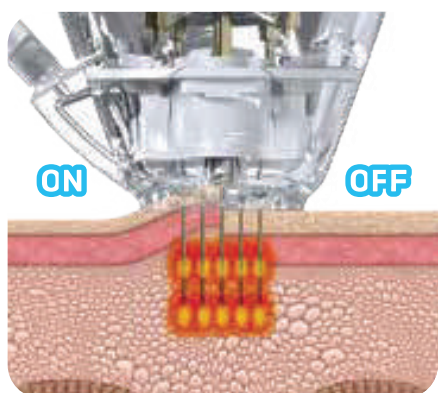

- ✓ Prevents skin folding or pushing during treatment
- ✓ Ensures precise needle-skin contact
- ✓ Guarantees accurate penetration and RF emission
- ✓ Activates microneedles only at optimal vacuum levels
- ✓ Easy ON/OFF switch

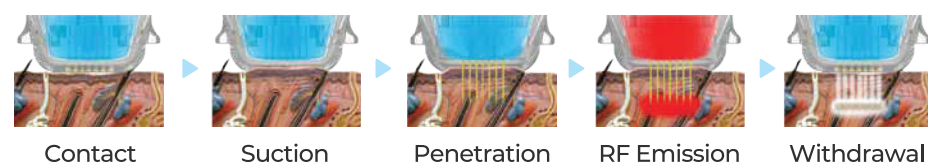

Dual Energy


1 MHz

Deeper Penetration, Enhanced Collagen Stimulation, Improvement in Skin Laxity, Longer-lasting Results. Beneficial for targeting loose and saggy skin, deep wrinkles, and lower eyelid fat

2 MHz

Surface-level Treatment, Precise Control, Quicker Recovery. Targets the upper layers of the skin, suitable for scars, enlarged pores, and acne.

Gold-plated microneedle

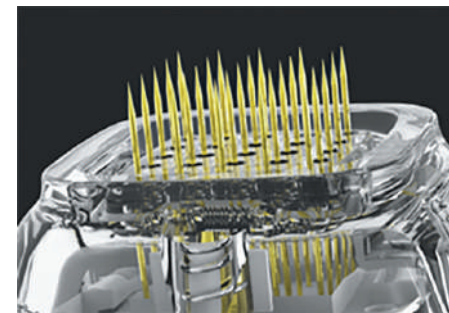

The gold-plated insulated microneedles minimize allergic reactions that could occur after treatment.

Auto parameter



Simply click on the area you want to perform, and the parameters for that area will automatically be displayed and help with easy and quick treatment.

MPC (M PURE CELL)



M


Human Bone Marrow-Derived Clonal Mesenchymal Stem Cell

PURE

High-Purity


CELL

Stem Cell Conditioned Media

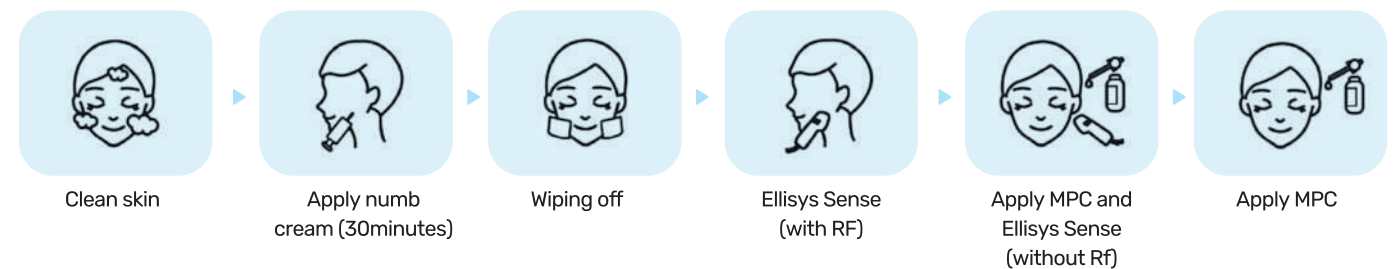

Various needle tips

Face

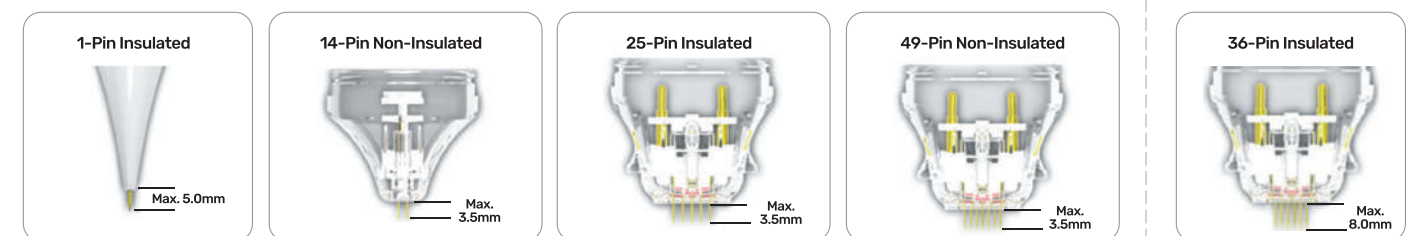

Body