




 Leaseir
AESTHETIC INTELLIGENCE



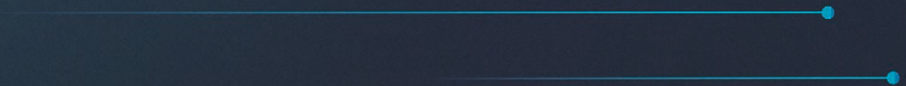


**We are experts
with over 20 years
at the forefront
of technology.**

Borned in Spain, where we have our facilities for production and R&D, we are committed to **quality** and **service**.



 Leaseir
AESTHETIC INTELLIGENCE

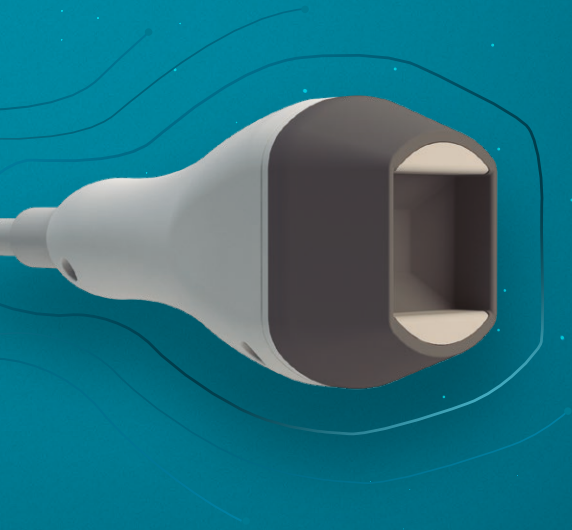

Leaseir is quality in treatments

LASER

Specializing in **diode laser** technology **for over 20 years**, our mission is to change the paradigm of hair removal worldwide, making this treatment **efficient, safe**, and **accessible** for everyone.



RADIOFREQUENCY

Our expertise in the beauty industry allows us to **innovate with facial and body devices** that continues to meet and exceed our clients' expectations.



Leaseir is innovation in products



MHR Diode Laser

The most powerful hair removal diode laser on the market.



DS Diode Laser

Portable diode laser for hair removal.



SRF Radiofrequency


Bipolar radiofrequency for facial and body treatments.



Leaseir MHR

the intelligent choice for hair removal specialists

THE **MOST
POWERFUL** DIODE
ON THE MARKET.

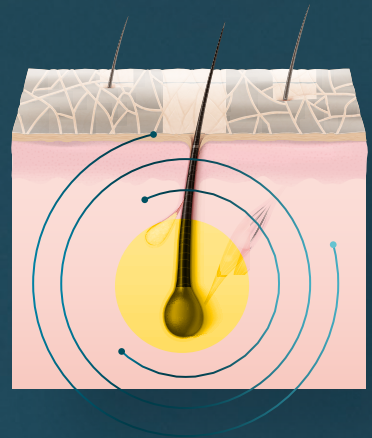

**REMOVES THE HAIR
FOLLICLE** WITHOUT
LEAVING RESIDUAL
HAIR.



+3000 DEVICES SOLD
since 2017



#1 hair removal **DEVICE** in Spain, Italy
and other countries.



Short pulses (from 4ms) and **high fluences** (up to 60J/cm²) to get permanent hair removal results in 4-6 sessions.



The **fastest device at maximum power** (4Hz and 3,7cm² spot in static mode).



Versatile for all skin types all year round (incorporates 3 wavelengths and SHR mode).



The **highest profitability on the market** (no consumables or maintenance).



Easy to use, safe and **comfortable** thanks to the Smart-cryo cold technology (leading edge at 0°C during treatment).



Made in Spain



Made in EU



+ **45 countries**

"Reach out to our experts to learn about the certifications of our devices."

HOW DOES IT WORK?

**We are the #1 choice
of specialists.
Our secret is the
combination of quality
and speed**



**High fluence
(up to 42J/cm²)**

To be able to remove the hair follicle without leaving residual hair.

**Short pulses
(from 4ms)**

Concentrate the energy in a few milliseconds and be more effective.

**Cold system
with 0°C tip**

Allows intensive but safe and comfortable treatment.

**HIGH
QUALITY**



**HIGH
SPEED**

**Spot area
(3.7 cm²)**

Covers more skin surface of skin with each shot.

**High frequency
(up to 4Hz and SHR)**

More shots per second at full power.


No breakdowns

Works 24/7 every day at maximum intensity with low downtime and low ratio of breakdowns.

Why is Leaseir MHR the device you need?

A **good treatment** has a **good balance** between short pulses and energy. It is about removing the hair follicle, NOT weakening the hair.

3 ways to give the same energy (pulse duration):

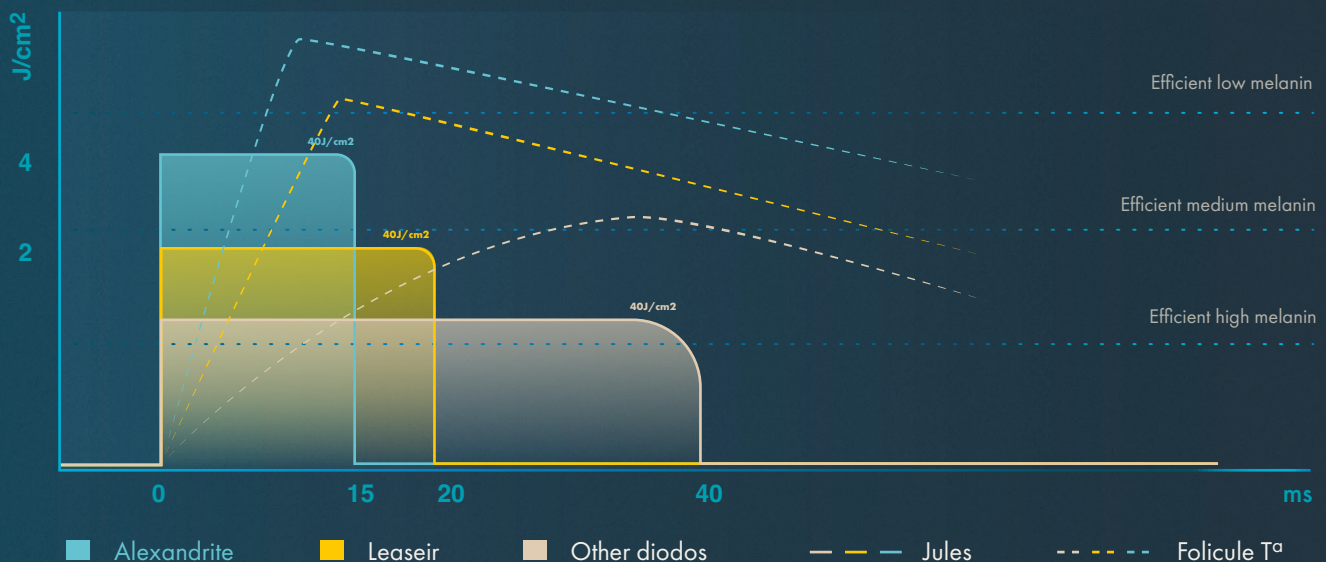

Cost effective device does **treatments in a short time** to allow you to do more treatments per day.


STATIC MODE			DYNAMIC MODE (SHR)		
4HZ (1Hz ankles and knees)			10HZ		
1 pass across the skin			4 passes over the skin		
Single Spot 1.17cm ²	Dual Spot 2.0cm ²	Quad Spot 3.7cm ²	Single Spot 1.17cm ²	Dual Spot 2.0cm ²	Quad Spot 3.7cm ²
4cm ² /sec	7cm ² /sec	14cm ² /sec	2cm ² /sec	4cm ² /sec	7cm ² /sec
11 minutes	8 minutes	5 minutes	16 minutes	12 minutes	9 minutes

1 female legs
session of 4.000cm²t

Our device combines both factors.


Comparison for a real hair removal session:

Fitzpatrick phototype 2-3 / Pulse 30ms / Maximum fluence / Maximum speed.

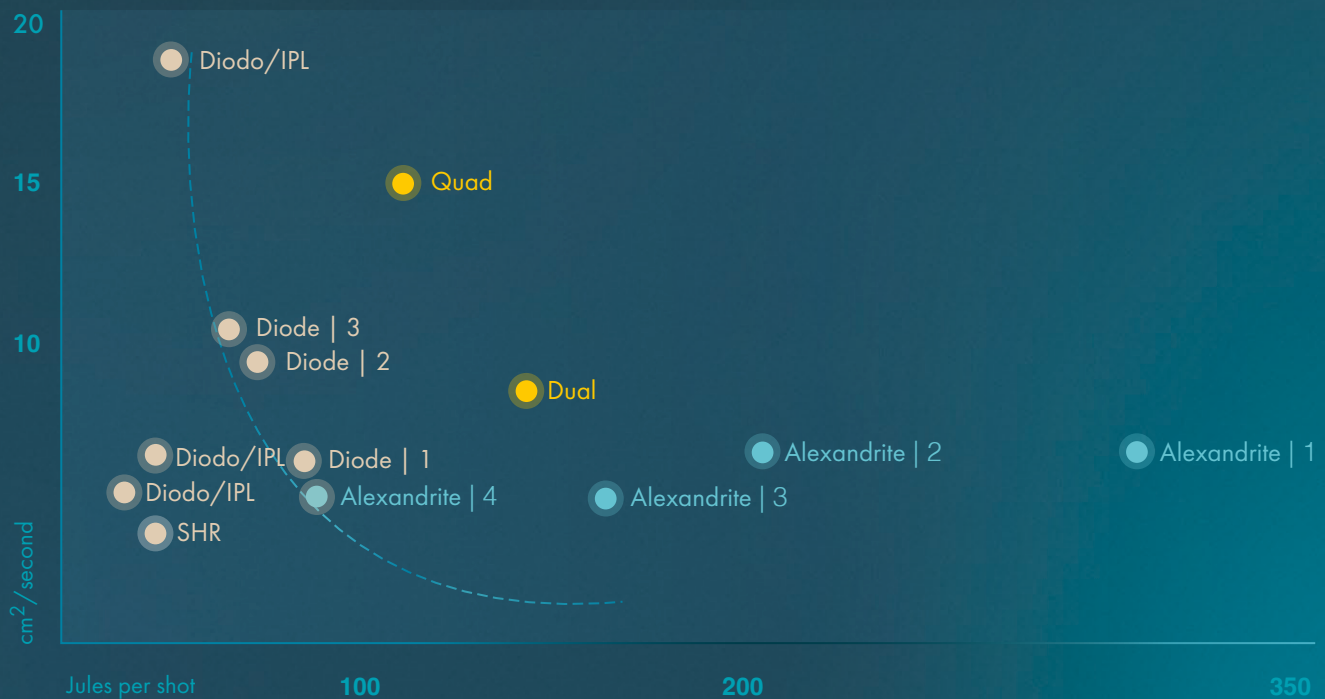

Designed for high profitability...

We simulate and compare the Return on Investment among different technologies, maximizing agendas at an average ticket of 50€.

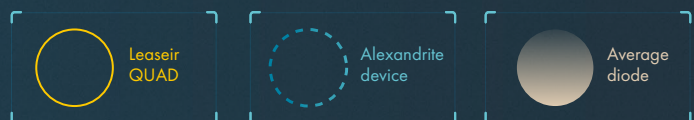

We become your partners to grow your business:

SERVICE 24/7

**LOW RATIO OF DOWNTIME
NO CONSUMABLES**


LEASEIR'S AI CLOUDMED
Online platform with your equipment real-time data.

...with high customer and user satisfaction


Before 1st session



Before 2nd session

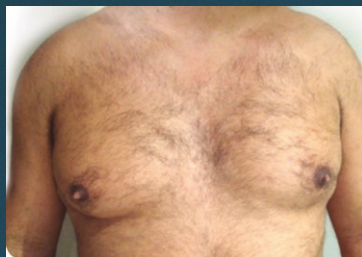

Before 3rd session

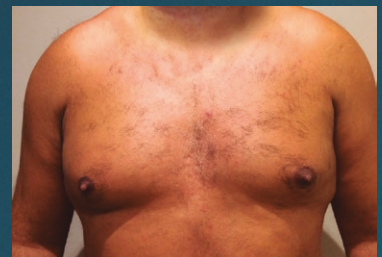

We acquired the Leaseir Single device at Dermoestética Affectum with excellent results. It stands out for its efficiency, safety, and a cooling system that minimizes discomfort, offering personalized and safe hair removal.

José Jiménez Manager
Dermoestética Affectum



We incorporated the Leaseir MHR device for its high speed, power, and effectiveness in hair removal, reducing hair density and making the process painless thanks to its cooling system.

Dr. Heisman Pereira. Manager
Clínica Pereira



We are the
best partner
for franchises
and specialized
chains.

SAFETY PROFITABILITY TRUST

Chain
#1 2 countries
450 MHR devices

Chain
#2 3 countries
240 MHR devices

Chain
#3 2 countries
145 MHR devices

Chain
#4 1 country
75 MHR devices

Characteristics



Console

Tip temperature	0-3°C
Dimensions	110cm x 45cm x 57cm
Weight	53 kg
Ambient conditions	15-25°C 0-90% Humidity 54-101.3 KPa
Laser classification	IV
Voltage	100-130V & 220-240V

Laser

SPOT	Single	Dual	Quad
Size	1.17 cm ²	2.0 cm ²	3.7 cm ²
Wavelengths	775 nm 810 nm 810+940+1060 1060	775 nm 810 nm 810+940+1060 1060	775 nm 810 nm 810+940+1060
Máx. Fluency	42J/cm ²	42J/cm ²	42J/cm ²
Pulse length	4-400ms	4-400ms	4-400ms
Frequency (static mode)	1-4 Hz	1-4 Hz	1-4 Hz
Frequency (SHR mode)	10 Hz	10 Hz	10 Hz
Consumables	No	No	No

About us

Leaseir is a company born in Spain, where we design and manufacture our own technology with passion and hard work. Our team has been at the forefront of diode laser technology for **more than 20 years**, making us one of the **world leaders**.

We are specialists and our goal is to change the paradigm of hair removal in the world, making this treatment **efficient**, **safe** and **accessible** to everyone.





HOW TO CONTACT LEASEIR:



Jimena Fernandez de la Vega 140 1G,
Gijón (Asturias). Spain



+34 984 39 09 34



info@leaseir.com